

PERMEADERM®


Biosynthetic Wound Matrix^{1,2}

See your patient through healing with flexibility and convenience with a proven next-generation technology^{1,2}





Easy to Use	Convenient and Versatile	Proven Effective
<ul style="list-style-type: none"> • Simple, one-time application • Transparency allows at-a-glance wound assessment • Quickly and intimately adheres to wound surface • Flexible and adherent even on articulating joints • Highly water resistant 	<ul style="list-style-type: none"> • Suitable for wounds, donor sites, and over meshed autografts • Expanded treatment utility with large sheets (up to 2,903 cm²) and glove options • Less prone to stick to secondary dressings • Usable without a secondary dressing after adherence 	<ul style="list-style-type: none"> • Creates optimal healing environment allowing epidermal migration • Customizable moisture management (defined porosity through stretch) • Provides bacterial barrier • Protects wounds from trauma

Unique, patented, dual layer, biosynthetic wound matrix aids in the function of lost epidermis until re-epithelialization^{1,2}

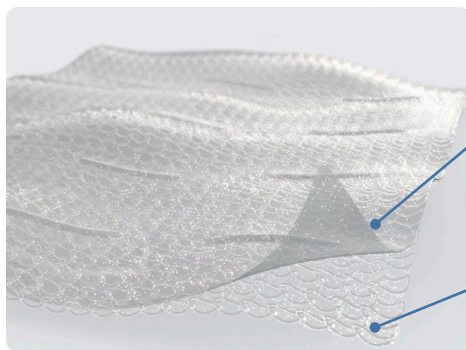

Fenestrated Outer 2D Silicone Epidermal Analogue facilitates diffusion of O₂, CO₂, and exudate to/from the wound, without sticking to a secondary dressing


Inner 3D Micro-Support System with biocoating, which expedites adherence to wound bed and supports early mobilization

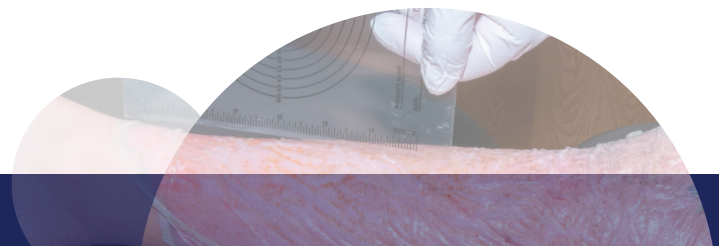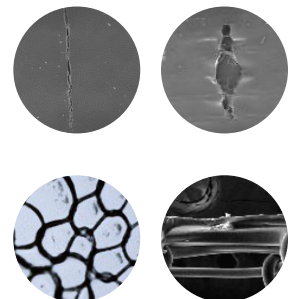

Case study: Deep Partial-thickness Burn (lower extremity)

12 year-old female • PermeaDerm® Sheet

Healing achieved within 15 days

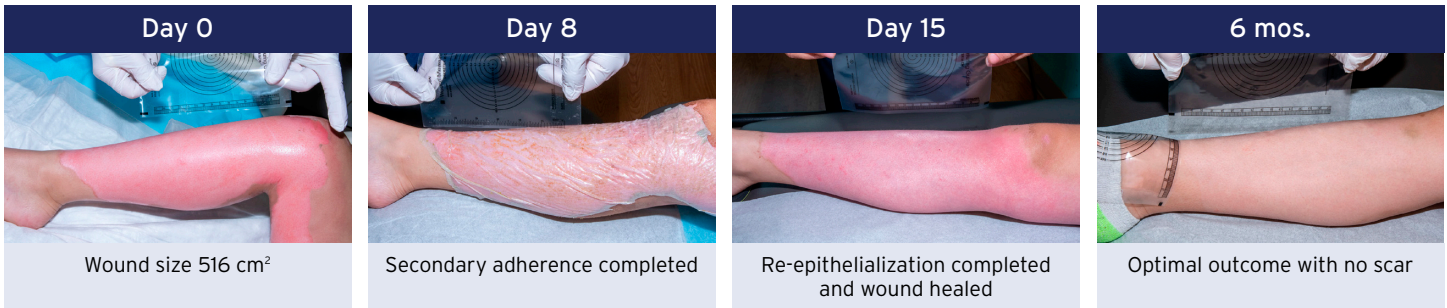

Case study courtesy of David Greenhalgh, MD FAC, Shriners' Hospital for Children.

Product Code	Product Description	Typical Use	cm ² /unit (per carton)	Units/ Carton
PermeaDerm® C				
C2050500	PermeaDerm 5" x 5"	Smaller wounds, such as chronic wounds, surgical wounds, trauma wounds, draining wounds	161 cm ² (1,610)	10
PermeaDerm® B				
B1051000	PermeaDerm 5" x 10"	Large wounds such as partial-thickness burns, trauma wounds, surgical wounds, donor sites, coverage over meshed autograft	323 cm ² (2,584)	8
B2101500	PermeaDerm 10" x 15"		968 cm ² (2,904)	3
B4153000	PermeaDerm 15" x 30"		2,903 cm ² (2,903)	1
PermeaDerm® Glove (unit = 1 glove, gloves are individually sealed)				
G1XS0000	PermeaDerm Glove XS	Debrided partial-thickness hand burns	406 cm ² (813)	2
G2S00000	PermeaDerm Glove S		562 cm ² (1,125)	2
G3M00000	PermeaDerm Glove M		713 cm ² (1,426)	2
G4L00000	PermeaDerm Glove L		801 cm ² (1,601)	2
G5XL0000	PermeaDerm Glove XL		947 cm ² (1,895)	2

Ask about ordering PermeaDerm® today.

Please contact your AVITA Medical representative or customer service at 833-GO-AVITA (833-462-8482)



PERMEADERM®

References: 1. Woodroof A, Phipps R, Woeller C, et al. Evolution of a biosynthetic temporary skin substitute: A preliminary study. *Eplasty*. 2015; Jul 20:15:e30.
2. FDA PermeaDerm 510(k) Premarket Notification.

PermeaDerm is a Stedical Scientific, Inc. product. AVITA Medical is the exclusive distributor of PermeaDerm in the United States of America.

©2024 AVITA Medical®. All rights reserved. PERMEADERM is a registered trademark owned by Stedical Scientific, Inc. PRO-1336-v1

